



SIMPLE Workshop 2017

Dominika Elmlund, Hans Elmlund

January 17, 2017

Contributors:

cyril.reboul@monash.edu

dominika.elmlund@monash.edu

hans.elmlund@monash.edu

Address:

Dept. Biochemistry and Molecular Biology

School of Biomedical Sciences

Monash University, Bldg. 77

Clayton, VIC, Australia, 3800

Webpage:

www.simplecryoem.com

1 What Will Be Covered

This part of the workshop will cover

1. DDD (Direct Detector Device) movie pre-processing
2. 2D analysis with PRIME2D
3. *ab initio* 3D reconstruction from class averages using PRIME3D

with the SIMPLE program package. We will use the latest unreleased version of SIMPLE, available to the workshop participants for downloading and beta-testing on the SIMPLE homepage <http://simplecryoem.com/beta.html>.

2 DDD Movie Pre-Processing

2.1 Motion Correction and Dose Weighting

During image acquisition, beam induced motion, charging and stage drift introduce blurring in the integrated movies. The individual frames need to be aligned with respect to one another in order to restore high-resolution information. Because of radiation damage, the electron dose is kept low during image acquisition, resulting in a very low signal-to-noise ratio (SNR) of the individual movie frames. Until recently, our understanding of damage processes in cryo-EM was based primarily on results obtained with 2D protein crystals (REFS 44,45 REVIEW). The use of 2D crystals provides direct means for assessing the degree of damage by looking at the rate by

which high-resolution diffraction spots are fading as a function of the accumulative electron dose. The disadvantage here is that it is impossible to separate spot fading due to loss of crystalline order from intrinsic molecular damage. Images of single rotavirus VP6 particles recorded with a total exposure of 100 electrons/Å² indicated that single-particles are more robust to damage at low and intermediate resolution than previously thought (REF 46 REVIEW). In the same study, a dose-weighting scheme was introduced that maximizes the SNR of the integrated movie in a dose-dependent manner. This dose-weighting strategy has been adopted by most software packages and proven to be a powerful addition to the ever-growing set of single-particle analysis tools.

SIMPLE implements a program called [unblur_movies](#) for simultaneous motion-correction and dose-weighting, with the objective of maximising the SNR of the integrated movie in a motion- and dose-dependent manner. Input is a text file `filetab`, simply listing the individual movies in `*.mrc` format, the pixel-to-pixel (or sampling) distance `smpd` (in Å), the number of CPU threads to use `nthr`, the body of the outputted files `fbody`, the exposure time of the movie `exp_time` (in s), the dose-rate at which the movie was acquired `dose_rate` (in e/Å²/s) and the acceleration voltage of the electron microscope `kv` (in kV), in addition to a number of other optional parameters that we will not be concerned with here. If we execute `simple_exec prg=unblur_movies` in the prompt, we obtain brief instructions for how to run the program:

```
@!#> simple_exec prg=unblur_movies
USAGE:
bash-3.2$ simple_exec prg=simple_program key1=val1 key2=val2 ...
```

REQUIRED

```
filetab = list of files(*.txt/*.asc)
smpd    = sampling distance, same as EMANs apix(in Å)
```

OPTIONAL

```
nthr      = nr of OpenMP threads{1}
fbody     = file body
lpstart   = start low-pass limit{15}
lpstop    = stop low-pass limit{8}
trs       = maximum halfwidth shift(in pixels)
exp_time  = exposure time(in s)
dose_rate = dose rate(in e/Å2/s)
kv        = acceleration voltage(in kV){300.}
pspecsz   = size of power spectrum(in pixels)
numlen    = length of number string
startit   = start iterating from here
scale     = image scale factor{1}
frameavg  = nr of frames to average{0}
tomo      = tomography mode(yes|no){no}
```

Output consists of four files per movie (`movie1.mrc` in this example):

1. `<fbody>_intg1.mrc` is the frame-weighted, motion-corrected and dose-weighted sum
2. `<fbody>_forctf1.mrc` is the un-weighted sum of the aligned individual frames
3. `<fbody>_pspec1.mrc` is the power-spectrum of the uncorrected unweighted movie sum (left) and the corrected weighted movie sum (right)
4. `<fbody>_thumb1.mrc` is a down-scaled version of `<fbody>_intg1.mrc`

The `*forctf*` output is created because, even though frame- and dose-weighting improves the SNR, it may degrade the contrast transfer function (CTF) signal and this file will therefore later

be used to determine the CTF parameters of the integrated movie. The **intg** image is the one you will use for particle picking and extract your identified individual particle images from. The **pspec** and **thumb** outputs are diagnostic. Please, check so that the Thon rings are concentric and have similar radial intensity distribution in any given resolution shell after alignment (right part of power-spectrum)—if not, trash the image. Please, execute:

```
@!#> cd ~/workshop/SIMPLE
@!#> ls
1_DDD-movie-preproc/ 2_PRIME2D/          2_PRIME3D/
```

You see that we have one directory for each step. Each folder contains a subfolder **data** that contains the experimental cryo-EM data that we will process in addition to its associated parameters.

```
@!#> cd 1_DDD-movie-proc/
@!#> ls
data/
@!#> ls data/
info.txt  movie1.mrc movie2.mrc
```

The `info.txt` file lists the parameters associated with the data.

```
@!#> cat ./data/info.txt
molecule: proteasome
exp_time=7.6s
dose_rate=7e/A2/s
kv=300
cs=2.7
smpd=5.26
frac=0.1
```

These movies are of proteasome molecules in ice, exposed for `exp_time=7.6 s` with a dose rate of `dose_rate=7 e/A2/s`. The electron microscope used to acquire these images has an acceleration voltage of `kv=300 kV` and a spherical aberration constant of `cs=2.7 mm`. We have downsampled these images for more rapid processing to a sampling distance of `smpd=5.26 Å` and we will assume 10% amplitude contrast when fitting the CTF `frac=0.1`. First, we need to create a file table, listing the movies to be processed with [unblur_movies](#).

```
@!#> ls data/movie* > movies.txt
```

Next, we execute the program.

```
@!#> simple_exec prg=unblur_movies filetab=movies.txt smpd=5.26 nthr=8
fbody=proteasome exp_time=7.6 dose_rate=7 kv=300
```

and we expect to see the following output in the terminal

```
>>> DONE PROCESSING PARAMETERS
>>> DONE BUILDING GENERAL TOOLBOX
>>> PROCESSING MOVIE:      1
>>> READING AND FOURIER TRANSFORMING FRAMES
100% |=====| done.
>>> WEIGHTED AVERAGE-BASED REFINEMENT
This % of frames improved their alignment: 100.
This % of frames improved their alignment: 100.
This % of frames improved their alignment: 100.
This % of frames improved their alignment: 100.
This % of frames improved their alignment: 47.
```

```

>>> LOW-PASS LIMIT UPDATED TO: 12.6667
This % of frames improved their alignment:    8.
>>> LOW-PASS LIMIT UPDATED TO: 10.3333
This % of frames improved their alignment:    8.
This % of frames improved their alignment:   53.
This % of frames improved their alignment:   42.
>>> LOW-PASS LIMIT UPDATED TO:  8.0000
This % of frames improved their alignment:    3.
This % of frames improved their alignment:    8.
>>> AVERAGE WEIGHT      :      0.0263
>>> SDEV OF WEIGHTS     :      0.0048
>>> MIN WEIGHT          :      0.0120
>>> MAX WEIGHT          :      0.0327
    50. percent of the movies processed

```

The first thing that happens is that all individual movie frames are read in and Fourier transformed. The algorithm is iterative and based on registration of the individual frames to a weighted average, where the weights are determined based on the agreement between the frames and the average. The algorithm automatically updates the resolution limit as the alignment accuracy improves, *i.e.* using only low-resolution Fourier components when the alignment errors are large and large movements need to be identified. It then successively updates the resolution limit as the alignment accuracy improves. When no improvements can be identified, the procedure stops. The average, standard deviation, minimum and maximum frame weight in the image series is reported upon completion of the alignment. The outputted image files are:

```

@!#> ls
data/                proteasome_intg1.mrc    proteasome_thumb1.mrc
movies.txt           proteasome_intg2.mrc    proteasome_thumb2.mrc
proteasome_forctf1.mrc proteasome_pspec1.mrc
proteasome_forctf2.mrc proteasome_pspec2.mrc

```

The first images we inspect are the power spectra. For visualisation, use `e2display.py`. Execute `e2display.py` in the current working directory, click on one of the `*pspec*` files and then Show 2D. To get an idea of the effect of the movie pre-processing in real-space, expand the

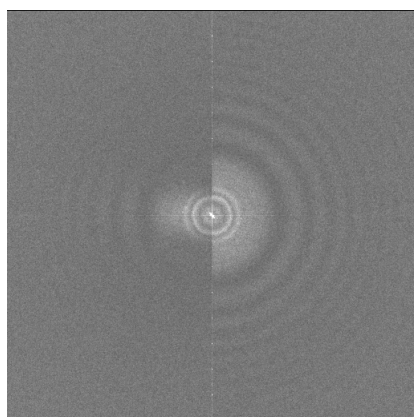


Figure 1: The left part of the power spectrum represents the movie integrated without applying the alignment parameters, whereas the right part is the aligned integrated movie. We see that the motion correction restores coherence, giving radially symmetric Thon rings extending to much higher resolution.

data folder from within `e2display.py`, then expand the `4comparison` folder and click on the `proteasome_straight_intg1.mrc` file, followed by Show 2D+. Do the same for the corrected movie `proteasome_intg1.mrc` and compare the two images.

2.2 CTF Parameter Determination

The contrast transfer function (CTF) mathematically describes how aberrations in a transmission electron microscope modify the image of the sample. The assumption we rely on is that the image

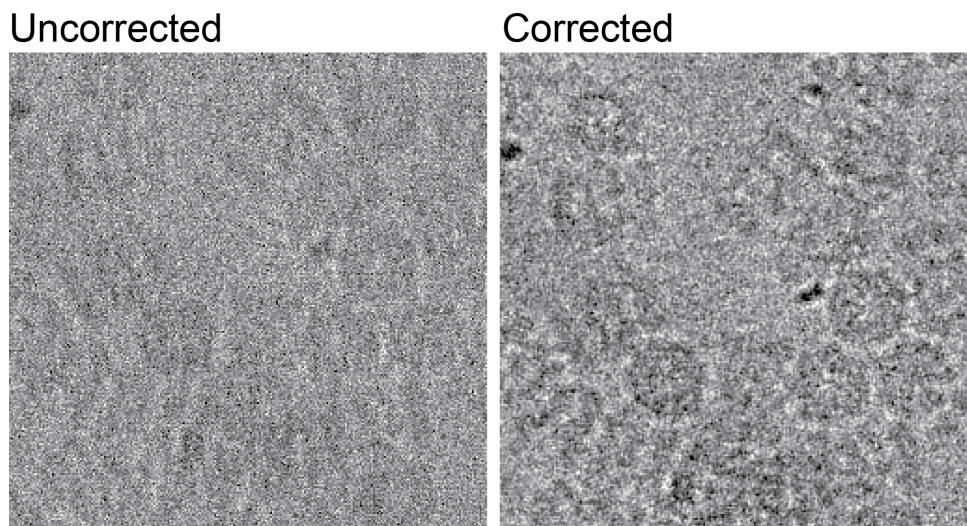


Figure 2: The uncorrected movie, integrated without alignment, is blurred and lacks high-resolution contrast, whereas the corrected movie has much higher contrast and better definition.

is a projection of the Coulomb potential distribution of the specimen, convoluted (multiplied in Fourier space) with the CTF and with Gaussian noise added. All 3D reconstruction methodology rely on this or variants of this assumption. The CTF is an oscillatory function (we saw the effect of it in the power spectrum in Figure 1) and the frequency of the oscillations depends on the spherical aberration constant cs in mm, which is an instrument specific constant, and the defocus in microns, which we systematically vary during image acquisition. The reason that we vary the defocus during data collection is that we want to avoid systematic lack of information in certain resolution regions, as no information other than noise is present in the zero crossings of the CTF. By varying the defocus we shift the positions of the zero crossings of the CTF. If astigmatism is present in the image, the Thon rings are no longer radially symmetric but elliptical and we need to take this into account as well.

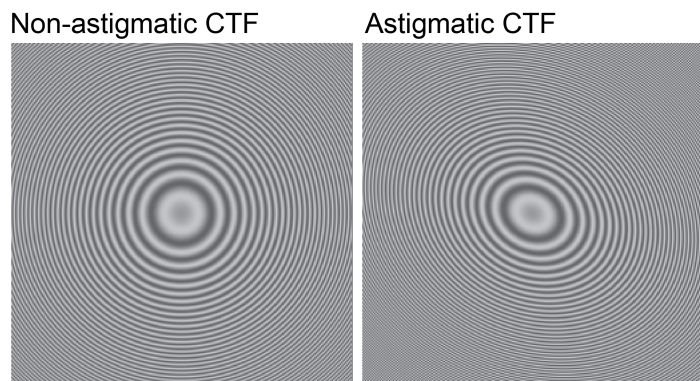


Figure 3: Non-astigmatic CTF (left) vs. astigmatic CTF (right)

SIMPLE implements a wrapper program called `ctffind` that executes CTFFIND4 (REF) version 4.1.X, producing a SIMPLE conforming CTF parameter document `deftab`. Input is a text file `filetab`, simply listing the individual integrated movies in `*.mrc` format.